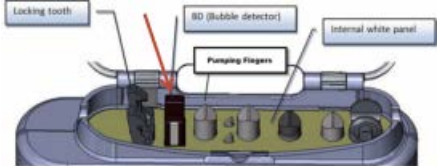


Alarms: When an alarm sounds, pump will immediately provide message with the error and suggested solution. Follow suggested alarm resolution, and then press OK key. To mute auditory alarm, press the **Mute** key.

Screen Display	Solution
Air in Line 	- Mute: Silences the auditory alarm for 2 minutes. - Check to be sure all line connections are tight. Check IV tubing for air bubbles. If there is air, disconnect IV tubing from patient and press the Prime key to remove the air from IV tubing, - Reconnect the IV tubing to the patient's catheter. - Clean Bubble Detector with a damp Q-Tip (use water only).  - Press OK key to resume infusion.
Insufficient Battery	Low battery voltage for current rate. Connect pump to power supply cord and wall outlet.
Battery Depleted	Pump will automatically shut down in 3 minutes. Connect pump to power supply cord and wall outlet.
Battery Reminder	End of battery life. Contact HomeMed to replace pump.

Alarms: When an alarm sounds, pump will immediately provide message with the error and suggested solution. Follow suggested alarm resolution, and then press OK key. To mute auditory alarm, press the **Mute** key.

Screen Display	Solution
Cassette Misplaced	The cassette on the IV tubing is not loaded correctly. Open the safety door, remove cassette, and reload the cassette.
Check for Occlusion	Verify all clamps are open and IV tubing is not kinked.
- Downstream Occlusion - Upstream Occlusion - Occlusion	- To clear occlusions check: - o All clamps are open; - o Cassette on IV tubing is inserted correctly; - o IV tubing is not kinked; - o No occlusion at the connection to your catheter. - o Once all occlusions are cleared press OK key to continue.
Flow Error	Check cassette on the IV tubing is inserted correctly. Check to be sure if battery is sufficiently charged. If not, connect pump to power supply cord and wall outlet. If alarm reoccurs contact HomeMed.
Internal Error	Pump will automatically shut down in 3 minutes. Contact HomeMed to replace pump.
Mechanism Error or Pump Fault	A pump fault has occurred. Contact HomeMed to replace pump.